



# WARMJET IV

AKAS's Whole new product line for  
the prevention and treatment of  
**Hypothermia**

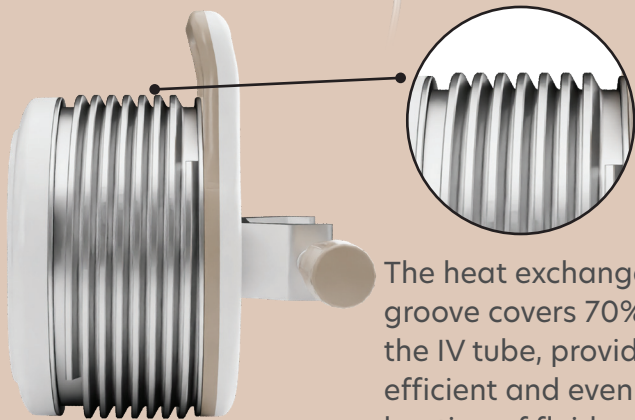

The WARMJET IV is designed to be the most compact and lightweight fluid warmer available, weighing less than 2 kg. You only need an IV set to operate the WARMJET IV. No extra disposables are necessary.

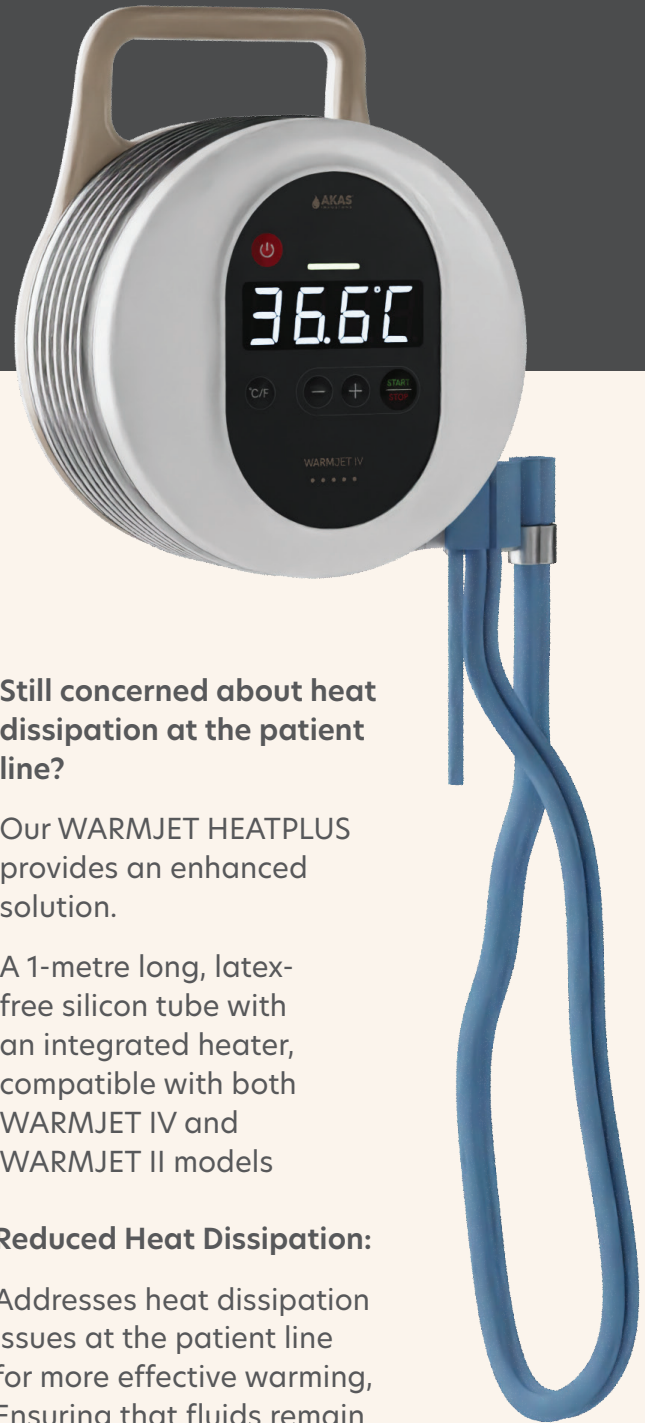

The heat exchanger groove covers 70% of the IV tube, providing efficient and even heating of fluids.

## WARMJET IV

### Specifications

Power Supply:	100-300VAC, 50/60Hz
Average Power Consumption:	150W
Maximum Power Consumption:	220W
Temperature levels:	36°C to 42°C
Temperature Increment:	0.1°C
Warm-Up Time:	2 min
Compatible IV line:	up to 5mm
Weight:	<2kg
Dimension:	183x162x210 mm
Flow Rate:	Up to 6000 ml/Hr

### ● Still concerned about heat dissipation at the patient line?

Our WARMJET HEATPLUS provides an enhanced solution.

A 1-metre long, latex-free silicon tube with an integrated heater, compatible with both WARMJET IV and WARMJET II models

### ● Reduced Heat Dissipation:

Addresses heat dissipation issues at the patient line for more effective warming, Ensuring that fluids remain at the desired temperature right up to the point of infusion.

### ● Dual Temperature Feedback:

Provides accurate and consistent warming with dual temperature feedback